

HEADLAMP

Illumination Systems

Light weight, bright, adjustable fiber optic headlight for medical and industrial applications

Adjustable glass fiber headlight with variable or fixed spot size. Comfortable, secure, adjustable high intensity light where you need it.



Rigid, adjustable headband and flexible crown adjustment make the headlamp assembly light weight and secure, or opt for the rigid, but fully adjustable model (not shown)

OPTIONS


- Single cable design
- 6', 7.5' or 10' cable length
- Totally flexible headband
- Integral input tip (customer specified)

FEATURES


- Adjustable, lightweight rigid head band with flexible (removable) or rigid (adjustable) crown support.
- Fixed 100mm spot or 10-75mm adjustable spot at 16". (406mm)
- Integral autoclavable joystick for accurate and easy positioning.
- Ergonomically designed lensing system is unobtrusive, yet perfectly positioned to put light wherever the user wants to see.
- 8 ft (2.43M) 5mm bifurcated multi strand glass fiber optic light guide for dependable long life.
- Threaded universal adapter accepts Storz, Olympus, Pilling, Wolf, ACMI and many other input tips. See reverse for full selection
- 1 year warranty.
- CE (IEC 601-1), UL(2601) certified as Class I equipment, type BF continuous operation.
- FDA approved to Class II specifications.

P/N	Weight	Cable Type	Description
FTIHL17524	.50 lbs (.22kg)	Bifurcated Glass	Rigid Headband/swivel headlamp/Variable Focus
FTIHL17525	.35 lbs (.16kg)	Bifurcated Glass	Rigid Headband/swivel headlamp/100mm Fixed Focus
FTIHL17526	.40 lbs (.18kg)	Bifurcated Glass	Semi Rigid Headband/swivel headlamp/Variable Focus
FTIHL17527	.25 lbs (.11kg)	Bifurcated Glass	Semi Rigid Headband/swivel headlamp/100mm Fixed Focus

Headlamp Features

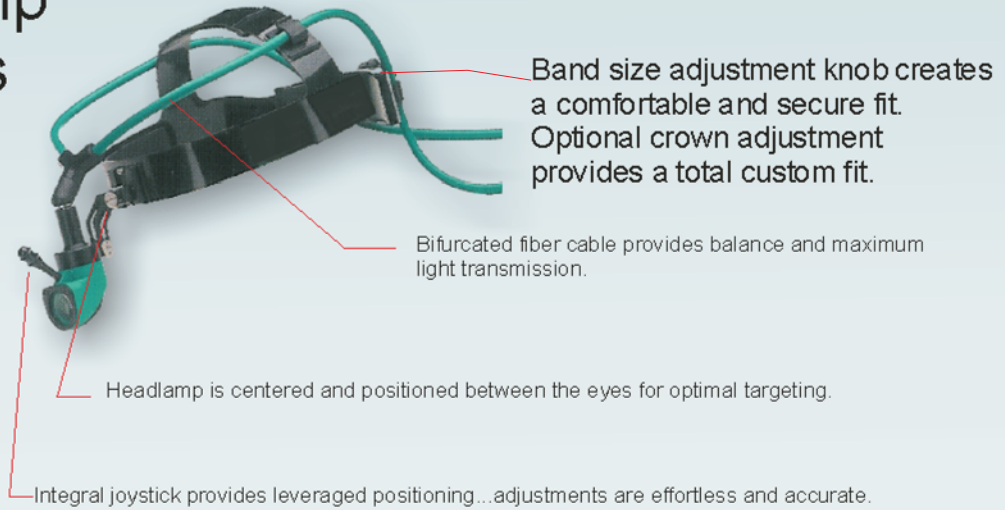


Input Fittings

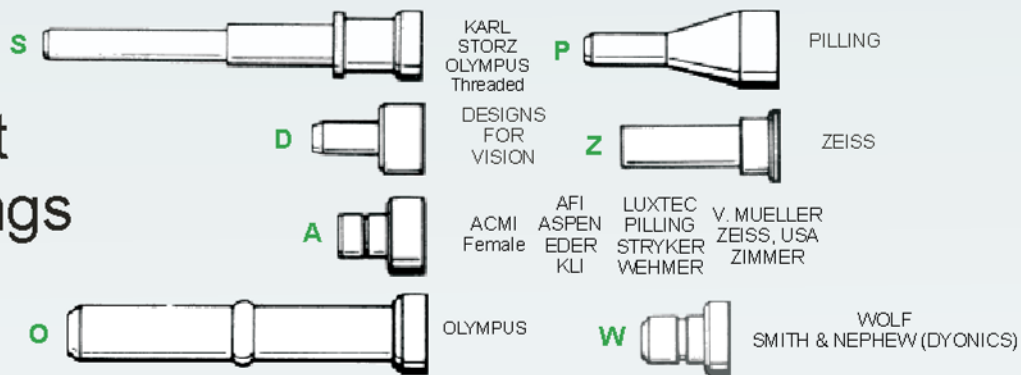

P/N	Letter code	Description
M13388	S	Storz chrome plated brass input
M13703	D	DFV chrome plated brass input
M13390	A	ACMI chrome plated brass input
M13389	W	Wolf/Dyonics chrome plated brass input
M13726	O	Olympus chrome plated brass input
M150002	Z	Ziess chrome plated brass input
SPECIAL	P	Pilling chrome plated brass input

The Pilling style endtip (Letter code P) cannot be made as a universal endtip. Requests for this style will be made as a custom item. Lead times apply.

Other types of end tips can be made to fit any light source. Call with requirements.

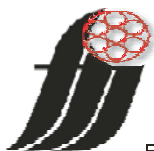

CORPORATE HEADQUARTERS
fiberoptics technology incorporated
ESTABLISHED 1977

1Fiber Road, Post Office Box 286, Pomfret, Connecticut 06258

PHONE 800.433.5248 • FAX 860.928.7884 • WEB SITE: www.fiberoptix.com • E-MAIL: info@fiberoptix.com